



GOLDEN



GOLDEN VERSIONS AVAILABLE:



CEILING

(code **GOLDENSO**)




MOBILE

(code **GOLDENPI**)



WALL

(code **GOLDENPA**)

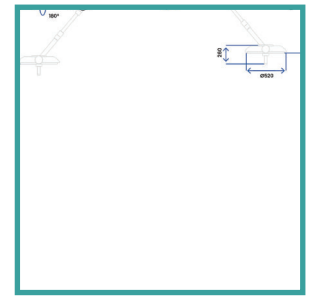

DOUBLE - CEILING

(code **GOLDENX2**)

GOLDEN

Golden is the first surgical lamp developed specifically for veterinary use. Usually, traditional surgical lamps are adapted for zoological use, but in this case Golden is the result of analysis and feedback from the veterinary sector.

Aware of the significant differences from the human sector regarding both design and use, the Rimsa R&D department has developed this innovative lamp. The Golden model is the only lamp studied, designed and produced for the veterinary needs.

WIDE LIGHT FIELD

The larger light diameter and the ability to use it at its best without the assistance of a medical team are just some of the features of Golden. The diameter of the light field can be adjusted between 11 and 35 cm giving the surgeon a very wide light field. A larger work area than a traditional surgical lamp allows the veterinary surgeon to focus during the operation instead of needing to move the reflector.

The adjustment of the diameter is possible through a capacitive sensor placed on the sterilizable handle, in order to guarantee the surgeon full control of the lamp.

ALWAYS ON FOCUS FUNCTION


The use of RIMSA technology with "always on focus" light guarantees perfect focus at any distance. In this way, the surgeon can freely move the lamp within the working area without having to worry about refocusing the light, which means having perfect light in every position.

HANDLING


The ease of movement allows to position the product directly by the surgeon, without the need of any assistant. This flexibility is possible thanks to the internal development of the arm structure and a perfect balancing of the reflector.




GOLDENX2

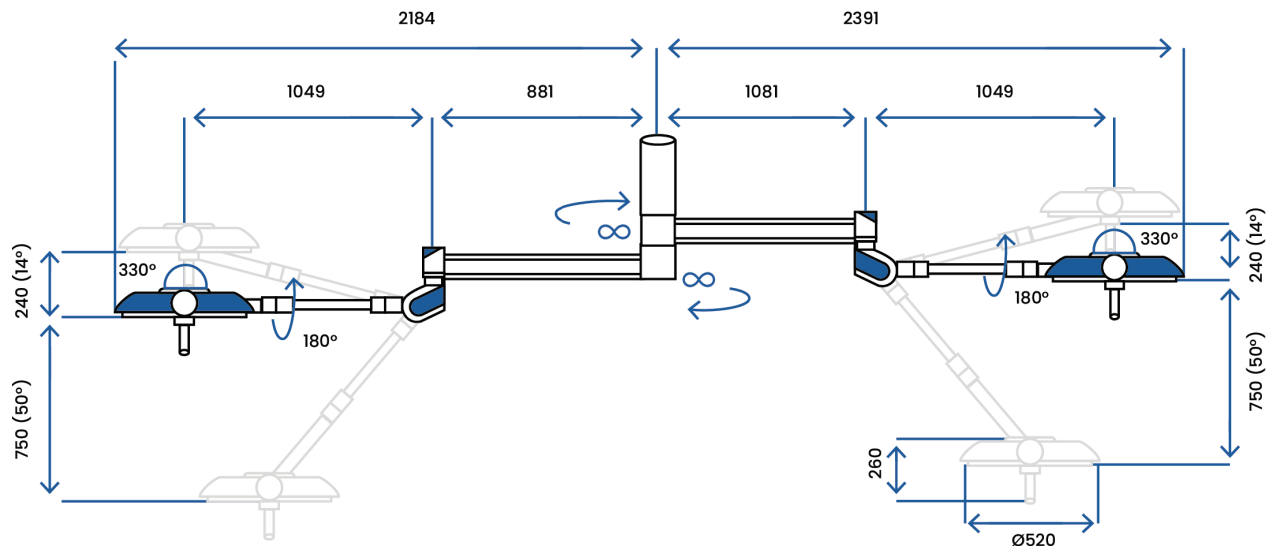


GOLDENPA

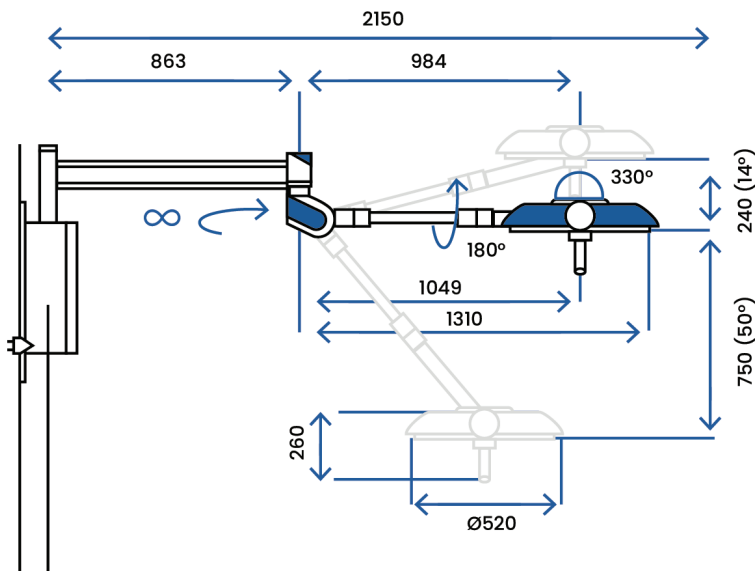

GOLDENPI

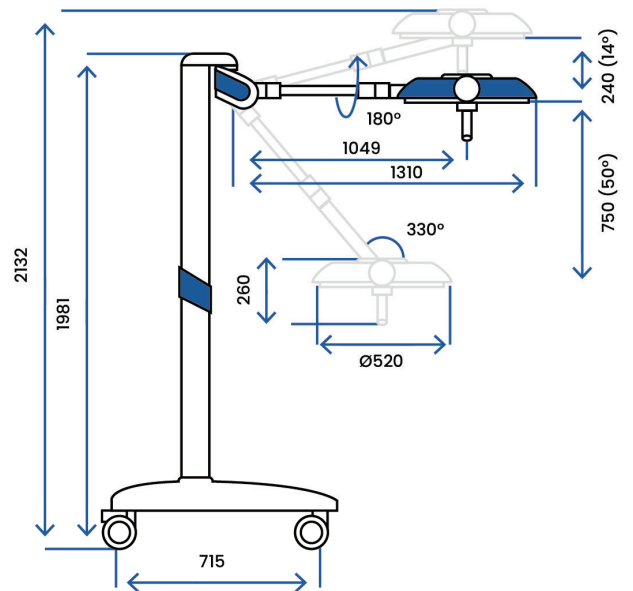

GOLDENSO

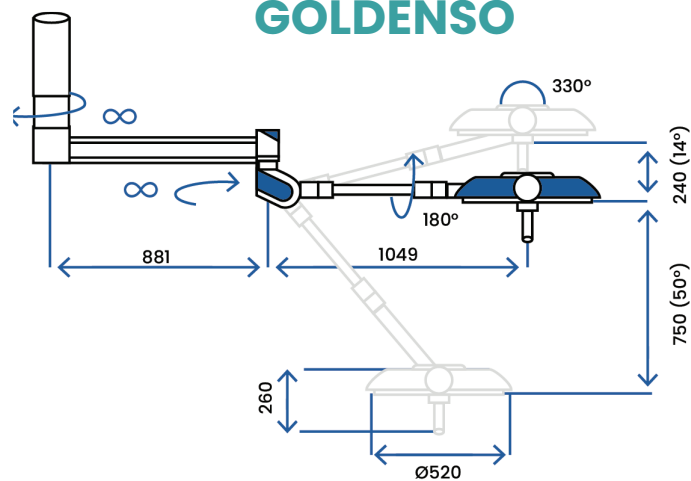

unit of measurement: mm



PERFORMANCES

Light intensity at 1 m distance (Ec)	Lux	160.000
Color temperature	K	4.000 - 4.600
Color rendering index (CRI)	Ra	94
Color rendering index R9		> 90
d10 light field diameter where illuminance reached 10% of Ec	mm	260
Diameter adjustment from - to -	mm	140-350
Depth of illumination IEC 60601-2-41 (L1+L2) at 60%	cm	81
Max irradiance	W/m ²	560
Irradiance / illumination	mW/m ² lx	3,47

POWER CONNECTION DETAILS

Primary alternate voltage (ac)	V	100-240
Frequency	Hz	50/60
Elettrical absorption	W - VA	60 - 67
Light head diameter	cm	52
N. of LEDs	LED	81
Average LED life	hours	> 60.000

